

HTP-1500 Heat Therapy Pump

Adroit
Medical Systems
800-267-6077

Thermoregulation & Temperature Therapy Products

HOSPITALS
for a
HEALTHY
ENVIRONMENT™ We Support H2e

adroitmedical.com

Patient side

- + effective pain relief
- + relax muscles
- + increase blood flow
- + promote healing
- + drug free

Bedside

- + tap water
- + mercury free
- + lifetime warranty*
- + real-time display
- + biotech friendly



HTP-1500 Heat Therapy Pump


Product Information

The HTP-1500 is the safest and most accurate localized heat therapy system available. Its bright and easy-to-read display provides continuous real-time temperature readings. The digital controller uses proprietary software to manage your prescribed set-point temperature to within one-degree Fahrenheit (1.0°F)!* Three temperature safety limits guard against the possibility of overheating and the Hi Limit Switches can be tested without taking the unit apart! The HTP-1500 is also Mercury Free! It uses Tap Water! And comes with a limited Lifetime Warranty!* Three standard sizes of heat therapy pads are available. Accessories include a rolling stand, bed bracket and temperature / flow rate tester.

Device Advise

Heat therapy is used to treat common pain disorders, including low back pain, pain from arthritis, trauma, muscle strain, wrist pain, and menstrual cramping. Continuous low level heat therapy (approx. 104°F / 40°C) is clinically proven to alleviate pain, relax muscles and increase tissue blood flow.¹ Both independent research and the FDA/CPSC recommend water-circulating systems over electric heating pads for consistent heating and better patient safety.²

¹ American Pain Society, Nov 2004
² FDA Public Health Advisory, 1996

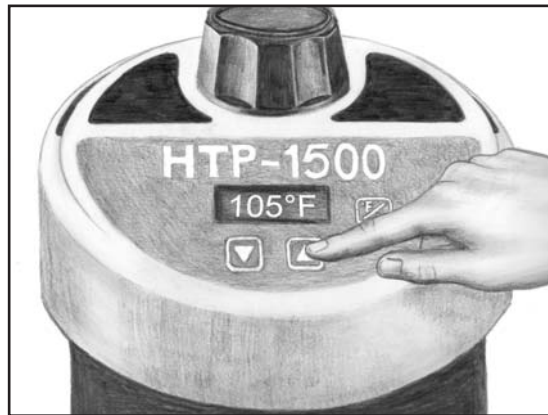

Easy to use touch-pad



Use the pad's blue side for dry heat or wet the white side for moist heat therapy**



Optional stand is perfect for bedside use
Order with or without a pad cradle

*1.0°F accuracy may not apply to units manufactured before 2005. Statement of warranty applies only to system components and workmanship and is not an implied warranty of fitness or suitability of the device for a particular purpose. Extended warranty beyond one year from date of purchase is discretionary and based on the hospital's continued purchase of heat therapy pads manufactured by Adroit Medical Systems, Inc.

**Never place a pad directly over an open wound without consent from the attending physician.

Pads & Accessories (more at web site)

ST-017	small pad 13.5" x 18.5" fabric both sides
ST-020	medium pad 16" x 22" fabric both sides
ST-024	large pad 20" x 26" fabric both sides
ST-220	home-care 16" x 22" fabric both sides
TFR-1	temperature / flow rate test unit
HTPS	stand with pad cradle
IVPB	I.V. Pole bracket
HTPB	bed bracket

Specifications

Set-point range.....	75°F(Off) to 107°F
Set-point accuracy*.....	1.0°F
Flow rate.....	≥15 gph
Reservoir.....	50 oz (1,500 ml)
Cord.....	18 AWG, 8 ft. Type SJT / hospital grade
Electrical.....	Volts 120, Hz 60, Amps 2, Watts 180
Approvals.....	UL 2601, CSA C22.2 No. 601-1 51DF